



Mubi Medical designs and manufactures healthcare solutions for the patient care environment. We focus on delivering thoughtful, practical solutions that are effective, budget conscious, and support patient recovery while maintaining high quality and aesthetics.

TO LEARN MORE ABOUT OUR SOLUTIONS, PARTNER WITH MUBI,
OR FIND YOUR AUTHORIZED DEALER VISIT:

www.mubimedical.com



**TRANSFORMING SPACES,
IMPROVING LIVES**

TS-700

S T R E T C H E R



The **TS-700** is a **VERSATILE** stretcher capable of supporting patient transport, examination, and recovery demands.



TRANSPORT



EXAMINATION
& RECOVERY



ER & TRAUMA



www.mubimedical.com

**Patient-ready,
700 lbs weight capacity**

Designed for Outstanding Mobility



Fold-away push handles with ergonomic grip to improve maneuverability and fold-away to provide access to the patient.

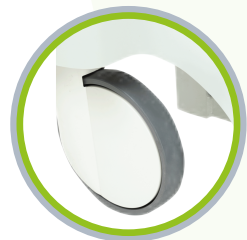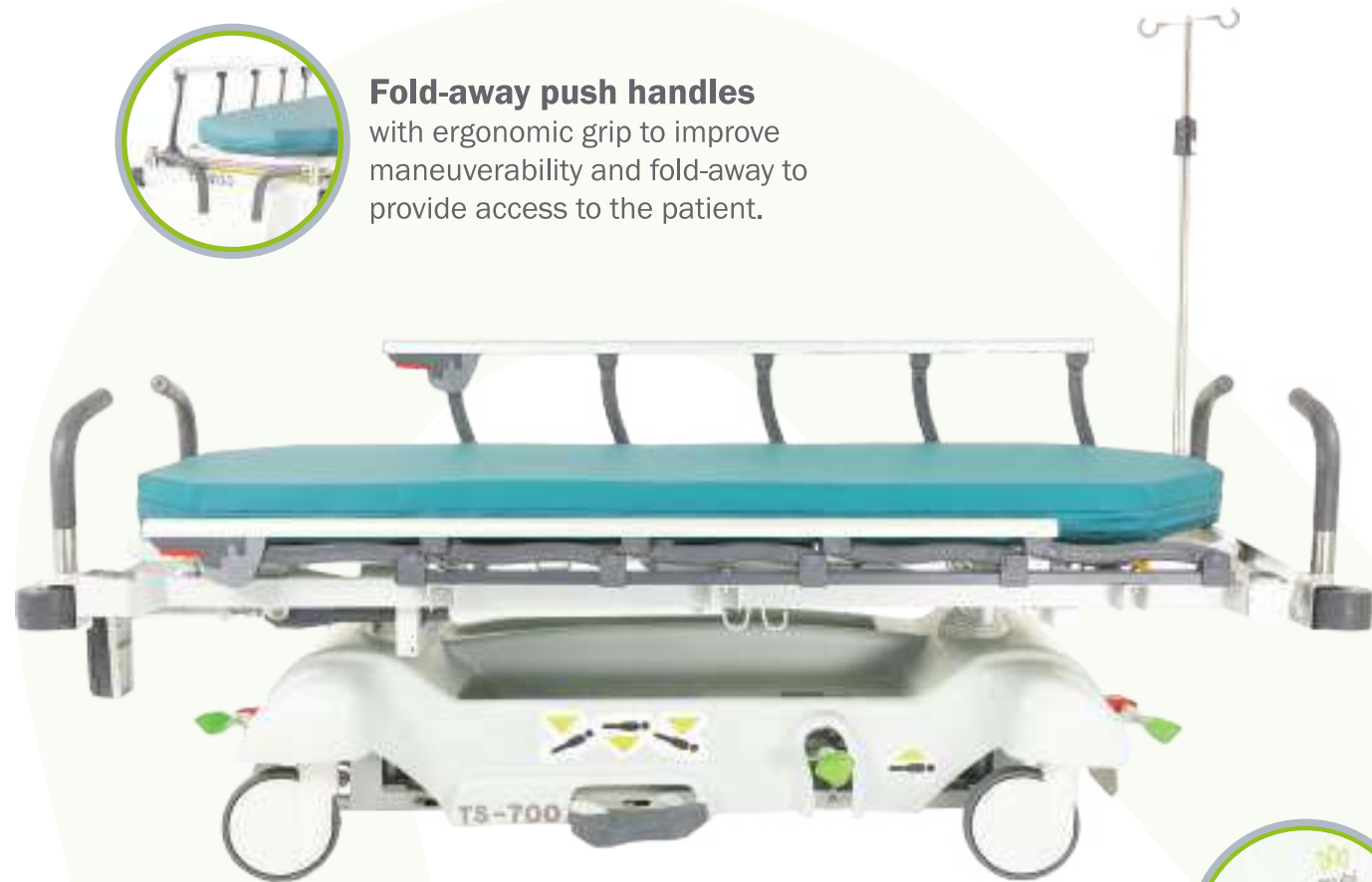

Dual-locking brakes prevent casters from rolling or swiveling as patients exit the stretcher.



Retractable fifth wheel allows optimal handling when making sharp turns, improving the stretcher control during transport.



Central locking breaks provide a stable platform for safe patient egress, improving fall risk prevention.

Designed for Safety

Recovery chair position

improves patient comfort, and facilitates treatments (airway clearance or administering medication).



Sliding Backrest reduces patient migration resulting in less boosting



Head of bed angle indicator aids with airway clearance head elevation compliance.



Range of patient redistribution surfaces with antibacterial and antistatic removable cover.



Knee gatch provides the patient added comfort. Easy to use manual crank.

700 lbs Safe Working Load

TS-700

STRETCHER

Designed for Quality and Durability

Heavy-gauge reinforced frame protected with rust resistant epoxy treatment and tested for 3000 hours in salt chamber.



Siderail bar, high strength aluminum, rust resistant and easy to clean.



Hydraulic elevation system protected with exceptionally durable and easy to clean anodized aluminum columns

Three-year Limited Product Warranty

Designed and Tested for 10 year expected life

TS-700

STRETCHER

Designed for Comfort and Ease of Use



Pneumatic assisted backrest reduces loads and helps prevent caregiver injuries.



Permanently mounted IV pole improves work flow by eliminating the hassle of missing IV poles.



Oxygen tank holder Integrated in base for storing O2 tanks sizes D and E.



Base provides ample room for storage of **patient belongings.**

Comfort Wide 28-inch deck

Feature-rich standard configuration offers outstanding mobility and safety.



Designed for Outstanding Mobility

- Push Handles, fold-away.
- Fifth Wheel, retractable.
- Central locking breaks.
- Dual-Lock Casters.

Designed for Quality and Durability

- Three-year limited warranty.
- Heavy-gauge reinforced frame.
- Protected hydraulic elevation system.
- Rust resistant siderails.

Designed for Comfort and Ease of Use

- Wider 28-inch deck.
- Integrated O2 tank holder.
- Patient belongings storage.
- Pneumatic assisted backrest.
- One-hand siderail release.
- Permanently mounted IV Pole.

Designed for Safety

- 700 lbs Weight Capacity.
- Recovery Position.
- Knee Gatch.
- Sliding Backrest.
- A range of Pressure Redistribution Surfaces.
- Head of bed and trend angle indicators.

TS-700

STRETCHER

TS-700

STRETCHER

Standard Configurations

- Included



	Transport	Examination & Recovery
2 section platform	●	
Fifth wheel	●	●
Fold-away push handles at head-end	●	●
Knee gatch / recovery position		●

Standard features (all configurations)

- 700 lbs weigh capacity
- 28-inch width patient deck
- Fifth wheel, retractable
- Fold-away push handles on head-end
- Sliding backrest
- 8-inch integral casters with wheel covers
- Pneumatic assisted backrest
- Trendelenburg and reverse trendelenburg
- Central-locking brake system
- Brake and steer pedals on head and foot ends
- Hydraulic elevation control pedals on both sides
- Fold down siderails with easy-to-clean aluminum tops
- 4 integrated IV sockets
- Drainage bag holders
- Bumper wheels on each corner
- Patient storage tray on base
- Integrated oxygen tank holder

Optional features

- Knee gatch
- Permanently mounted IV pole at head-end with 2 hooks
- Removable IV Pole, adjustable height and 2 hooks
- Fold-away push handles on foot-end
- 3-inch surface pad with antibacterial and anti-static cover
- 4-inch surface pad with antibacterial and anti-static cover
- Side Rail Pads
- Angle indicators for head of bed and trend
- Passive Leg Rest
- Monitor tray board
- Chart holder

Specifications

Weight Capacity 700 lb (318 Kg)
 Overall length 84" (215 cm)
 Overall width 34" (87 cm)
 Patient surface length 77" (195 cm)
 Patient surface width 28" (70 cm)
 Knee gatch 0° - 45°
 Backrest 0° - 90°
 Height range 23.25" (59 cm) - 34" (86.4 cm)
 Trendelengurg 17°
 Reverse trendelenburg -17°
 Caster diameter 8" (20 cm)

Warranty

Three-year limited product warranty

MUBI Medical reserves the right to change specifications without notice.