

OVERBED TABLE



PNEUMATIC
Ref. 84EU28



@mubimedical
www.mubimedical.com

COMFORT AND EASE OF USE

Its low profile of 1.96 in (50 mm) allows use under stretchers and beds with a 5th wheel, providing comfort and facilitating the work of healthcare staff.

DINING TABLE SURFACE

Made of durable material, free of sharp edges, and easy to disinfect.

Seamless design prevents fluid and contaminant buildup, allowing food and beverages to be placed safely without the risk of spills (up to 42 oz).

Stain, scratch, and disinfection product-resistant surface.



HEIGHT ADJUSTMENT HANDLE

Smooth one-handed ascent. To lower, press the release tab and apply gentle pressure to the surface.


CUP HOLDERS

2, with smooth edges that facilitate disinfection.



WHEELS

2 wheels with brakes on the column side.
2 wheels at the front.

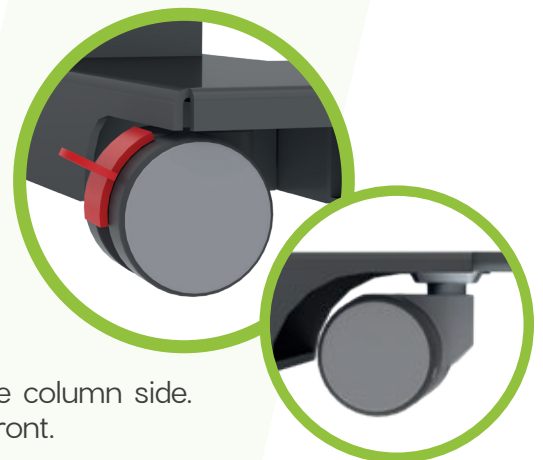

OPTIONAL ACCESSORIES

*With additional cost.



DRAWER

Located under the table, dual-access, for personal belongings, and with a flip-up mirror.

8.75 in L x 5.25 in W x 1.75 in H

Mirror: 9.1 in L x 5.6 in W



AUXILIARY BASKET


Provides storage for personal belongings while keeping them visible. Designed to be mounted on the side of the column.

1: Internal dimension (WxLxH): 5.3 x 9.3 x 9 in (134 x 236 x 228 mm)

2: Internal dimension (WxLxH): 4.3 x 4.3 x 4.7 in (109 x 109 x 119 mm)

SURFACE COLOR OPTION*

Surface table is available in a variety of wood grains and some solids.

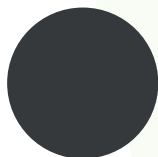

BLACK
PREMIERE
MATT



HALIFAX
OAK ABAK
SYNCHRO



HAYA
NATURAL



OAK MELINGA
GREY



PACIFIC B



WENGUE
REALISTIC
PORE



WINCHESTER
PA

* Color changes may cause changes in delivery time.

* Minimum production order of 10 units.

* Consult the guide and manufacturing times with your representative.

APPLICATION ENVIRONMENTS

Ideal for environments such as: recovery, observation, hospitalization, intensive care, and special care.

STANDARD FEATURES *

- Maximum Static Load Capacity: 100 lbs (45.3 kg), evenly distributed across the surface.
- Safe Functional Load Capacity: 55 lbs (25 Kg).
- High-precision thermolaminated tabletop with integrated cup holder.
- Slim-profile base designed for stability and compatibility with stretchers and beds equipped with a 5th wheel, without compromising maneuverability.
- Two dual-wheel casters with brakes, made of polyurethane, diameter: 2 in (50.8 mm).
- Two free dual-wheel casters, made of polyurethane, diameter: 1.4 in (35 mm).
- Assembly-friendly design optimized for fewer steps and reduced setup time.

MANUFACTURING AND FINISHING *

- Structure made of cold-rolled steel sheet and aluminum profile column.
- Finished with electrostatically applied, rust-resistant epoxy powder coating.

DIMENSIONS **

Base Length	32 in (812 mm)
Base Width	19 in (482 mm)
Minimum Height	28 in (711 mm)
Maximum Height	43 in (192 mm)
Surface (length x width x height)	32 in * 17 in * 0.70 - 0.80 in (812 mm * 432 mm * 18 - 20 mm)

*Units expressed in inches (in) and millimeters (mm). Tolerances: ± 0.5 in (12.7 mm)

PACKAGING AND PACKING

- Individual packaging option. Box reference 15216.
- Container: 20 STD 140 unds.
- Product weight: 40 lbs / 20 Kg.

NOTES

** MUBI Medical reserves the right to modify specifications without prior notice.

*** Color change can have an additional cost and may affect time of delivery; please consult with your representative.



TRANSFORMING SPACES, IMPROVING LIVES



Medellín – Colombia
Calle 2 # 50-185
(57) 604 285 2838

Cali – Colombia
Calle 18 # 118-85
(57) 602 489 6999

