

TS-700

TRANSPORT AND EMERGENCY STRETCHER



COMFORT AND DYNAMISM



@mubimedical
www.mubimedical.com

TS-700 protects your patients while protecting your investment, thanks to a transport, examination and recovery stretcher that meets regulatory standards, ensuring quality care.

TS-700

TRANSPORT AND EMERGENCY STRETCHER

DESIGN AND SAFETY

High resistance materials and easy disinfection finishes, compatible with products used in hospitals such as: hydrogen peroxide, quaternary ammonium, sodium hypochlorite and antiseptic alcohol**.

Sliding backrest*

reduces patient migration resulting in less boosting, eases burden on caregiver's back and helps reduce injuries made of high-strength steel, easy to clean and disinfect.


* Four-section platform (Ref.44EU2610):
one fixed and three movable.

* Two-section platform (Ref.44EU2611):
one fixed and one movable.



Corner bumper

allows protect the structure and building walls institutions from possible impacts.

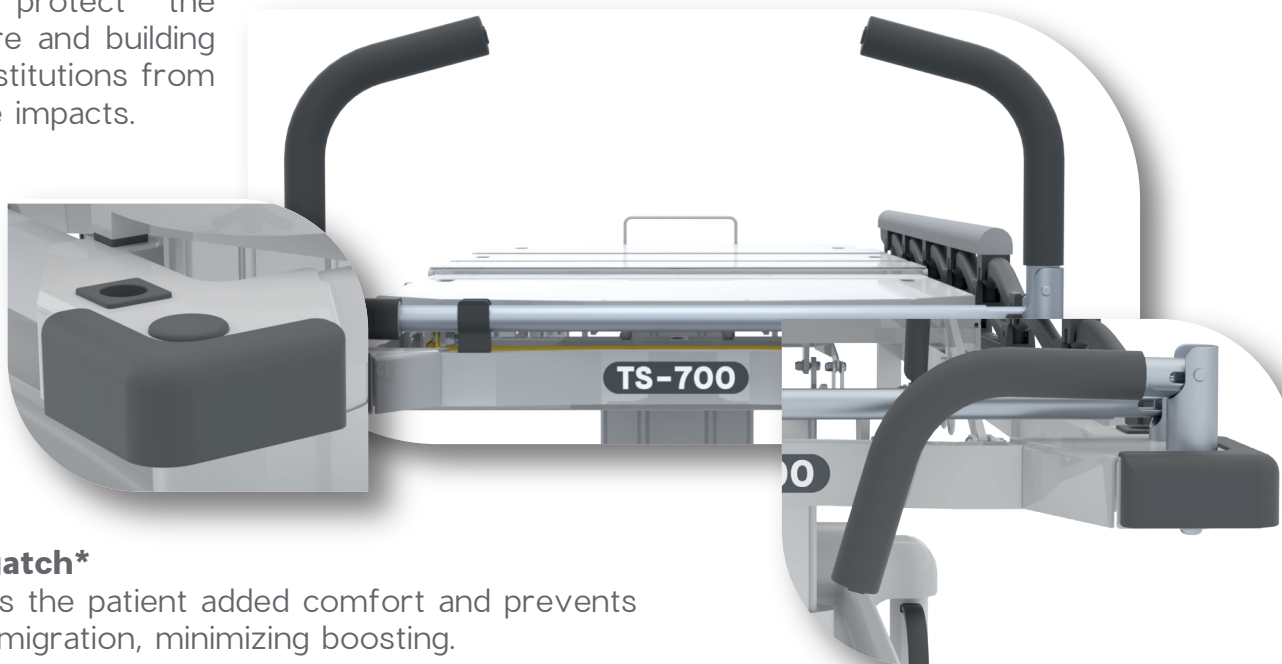

Knee gatch*

provides the patient added comfort and prevents patient migration, minimizing boosting.

* Only for 4-sections stretcher. Ref. 44EU2610



Fold-away push handles

facilitate transport, improve maneuverability and provide access to the patient.

** Recommended composition:

- Hydrogen peroxide 3% and 7.5% 200 ml peroxide per liter of water.
- Quaternary ammonium 4 ml per liter of water.
- Sodium hypochlorite 4 spoonfuls for each liter of water approximately.
- Antiseptic alcohol 60-80% 7 parts of alcohol x 3 parts of water.


IDEAL FOR EMERGENCY AND TRAUMA UNITS

TS-700

TRANSPORT AND EMERGENCY STRETCHER

SAFETY SIDE RAILS

One-section, foldable, made of high-strength polymer and aluminum cover, easy to clean and disinfect, simple one-handed operation, facilitating medical care.

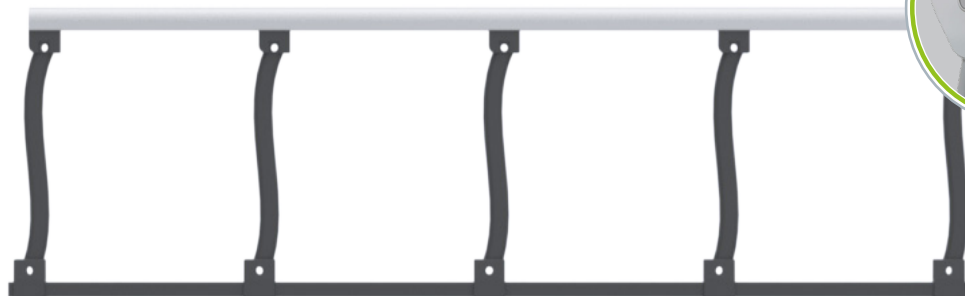

One-handed siderail release simple and efficient locking mechanism which avoids pinch point risks.



Central-locking brake system

prevent wheels from rolling or swiveling as safe patient ingress-exit.



Retractable fifth wheel

optimal handling, smooth and safe turns, greater control during transport, in blocked mode the wheel does not come into contact with the ground, allowing the infeed table to enter.


Improved fifth wheel adherence force against the ground, we upgraded the mechanism achieving greater force and adhesion.

TELESCOPIC IV POLE

One fixed and folding (left of backrest) and one removable to be placed in any of the 4 corners.



Removable IV pole
2 stage and 2 hooks.

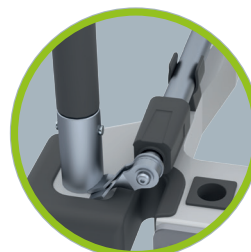

Permanently mounted IV pole
at head-end, 2 hooks.

TS-700

TRANSPORT AND EMERGENCY STRETCHER

*Ref. 44EU2611.



Head of bed angle indicator
aids with airway clearance head elevation compliance.



Hydraulic elevation system
exceptionally durable and easy-to-clean anodised aluminium columns protect the mechanism.



Base cover
for use as an integrated oxygen tank holder or as a space for patient belongings.

STANDARD DESIGN WITH MULTIPLE CONFIGURATION OPTIONS, WITH MOBILITY AND SAFETY IN MIND

TS-700

TRANSPORT AND EMERGENCY STRETCHER



COMFORT AND SAFETY

Ideal for hospital environments, high strength steel frame, foam and laminate surfaces; a precise combination of materials for support, comfort and reliability.

*Ref. 44EU2610.

SURFACE PAD (Optional with additional charge)

Helps in the treatment and prevention of pressure injuries, high density polyurethane (PU) foam core, $40 \pm 2 \text{ Kg/m}^3$ (3107.3 lb/ft³)*.

* foam thickness options (Optional with additional charge)

Ref. 740-533. 3-inch surface pad with antibacterial and anti-static cover.

Ref. 740-540. 4-inch surface pad with antibacterial and anti-static cover.



PAD COVER COLOR

TURQUOISE
HOSPITAL EXTREME


HOSPITAL EXTREME COVER

Made of nylon, permablok lacquer coating, light protection AATCC TM 16-3 / NTC 1479 / CFFA-2 a1 - Blue Wool, waterproof, anti-static ASTM D-257, bacteria resistance AATCC 147, JIS Z 2801, mildew resistant NO GROWTH (NG), invisible waterproofed seams, stitch resistant DIN 54301, flame retardant FMVSS 302 and latex free.

OVERALL DIMENSIONS *

TS-700

TRANSPORT AND EMERGENCY STRETCHER

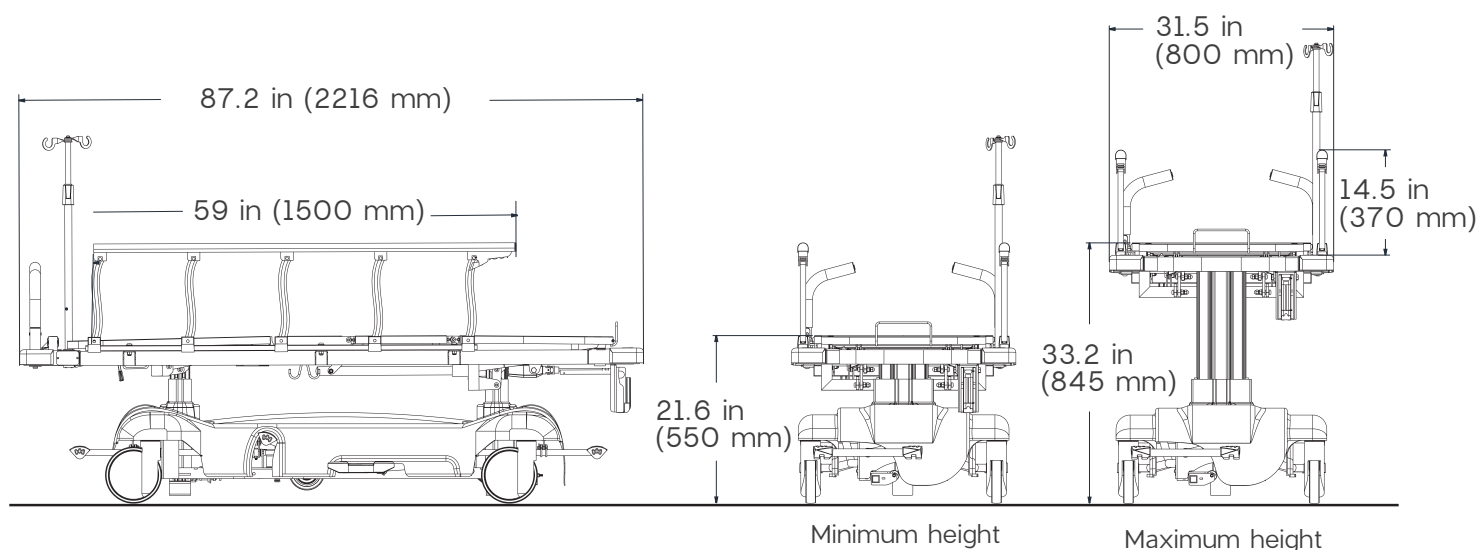

PATIENT SURFACE *

Length	77 in (1950 mm)
Width	26 in (660 mm)

MOTIONS * **

Backrest	0° a 89°
Footrest***	0° a 45°
Trendelenburg	15°
Reverse trendelenburg	-14°

* Tolerances: degrees $\pm 3^\circ$ dimensions 0.5 in (12.7 mm).

**Motions performed by pneumatic and hydraulic systems.

*** Only 4-sections stretcher. Ref. 44EU2610.

OPTIONAL ACCESSORIES (ADDITIONAL COST)

- Ref. 740-459-02. Fold-away push handles on foot-end.
- Ref. 740-533. 3-inch surface pad with antibacterial and anti-static cover.
- Ref. 740-540. 4-inch surface pad with antibacterial and anti-static cover.
- Ref. 740-461. Removable IV pole, 2 stage and 2 hooks.

MANUFACTURING AND FINISHING

- Structure made of high-strength cold rolled steel tubing, reinforced with structural tubing.
- Finished with electrostatically applied, rust-resistant epoxy-polyester powder paint.

PACKAGING

- Container: 40HC – 39 units, 40STD – 26 units, 20STD – 12 units.

TS-700

TRANSPORT AND EMERGENCY STRETCHER



*Ref. 44EU2611.

EXAMINATION AND RECOVERY CONFIGURATION



- Stationary foot patient platform
Ref. 44EU2611.

STANDARD FEATURES

- Retractable fifth wheel.
- Knee gatch / recovery position* (Optional with additional charge).
- Fold-away push handles on head-end.
- Sliding backrest.
- Permanently mounted IV pole at head-end, 2 hooks.
- Angle indicators for head of bed and trend.
- 8-inch integral casters with wheel covers.
- Pneumatic assisted backrest.
- Trendelenburg and reverse trendelenburg.
- Central-locking brake system.
- Brake and steer pedals on head and foot ends.
- Hydraulic elevation control pedals on both sides.
- Fold-down siderails.
- 4 integrated IV pole sockets.
- Drainage bag holders.
- Bumper on each corner.
- Patient storage tray on base.
- Oxygen tank holder on base.

TRANSPORT CONFIGURATION



- Includes knee gatch / recovery position
Ref. 44EU2610.

SPECIFICATIONS

- Weight capacity 700 lbs.
- Patient surface width 26 in.
- Patient surface length 77 in.
- Knee gatch 0°-45°*.
- Backrest 0° - 89°.
- Height range 21.6 in - 33.2 in.
- Trendelenburg 15°.
- Reverse trendelenburg -14°.
- Wheels diameter 8 in.
- Overall length 87.2 in.
- Overall width 31.5 in.

ACCESSORIES (additional charge)

- Fold-away push handles on foot-end.
- Removable IV pole, 2 stage and 2 hooks.
- 3-inch surface pad with antibacterial and anti-static cover.
- 4-inch surface pad with antibacterial and anti-static cover.

* Only for 4-sections stretcher Ref. 44EU2610

NOTES

- ACCESSORIES, please consult with your executive.
- MUBI Medical reserves the right to change specifications without prior notice.
- Limited Product Warranty, please consult with your executive.

TRANSFORMING SPACES, IMPROVING LIVES



Medellín – Colombia
Calle 2 # 50-185
(57) 604 285 2838

Cali – Colombia
Calle 18 # 118-85
(57) 602 489 6999



mubi

M E D I C A L

UPDATE DATE: 12/05/2025