



## ADVANCED CARE AND PERFORMANCE

**MaxCare**, your strategic ally — optimizing clinical workflows. It empowers caregivers to deliver high-quality care with greater efficiency and safety.

### PATIENT SURFACE WIDTH 26 INCHES

Stretcher with a two-section surface, made of high-resistance, easy-to-disinfect steel sheet.

Support for locating fastening straps in the back and feet section.

Mattress retainer on feet.



### ADJUSTABLE HEIGHT, TRENDLENBURG & REVERSE TRENDLENBURG

Ergonomic pedal system provides smooth control over height and Trendelenburg positions — improving caregiver efficiency and ease of use.


### BACK ADJUSTMENT

Reduces patient movement and the effort needed to reposition — enhancing patient comfort and protecting caregivers from strain and injury.




## SAFETY SIDE RAILS

Features 'tuck away' system, designed to minimize transfer spaces and protect patients during transport.

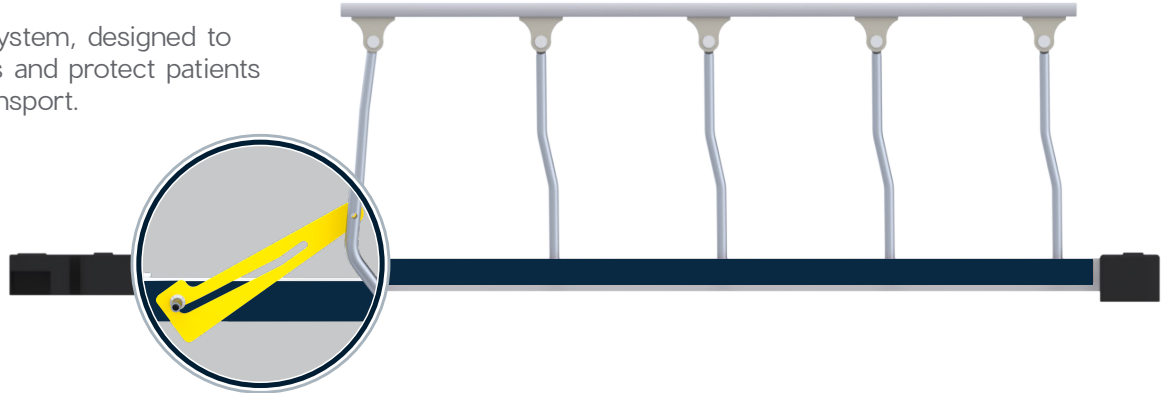

Easy-to-use release lock enables smooth folding of the collapsible side rail, improving patient safety, comfort, and caregiver efficiency

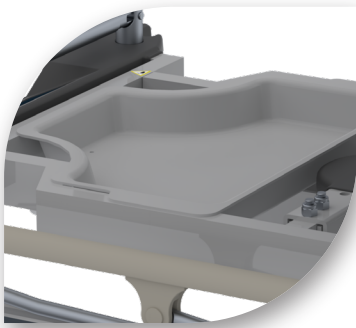

## AUXILIARY TRAY

Allows placement of medical supplies.



## FOLD DOWN PUSH HANDLES

Located in the back section, the push handles improve maneuverability and reduce caregiver strain during transport.



## BUMPER CORNER

Located in the corners to protect both the stretcher and walls from impact damage.

## PATIENT'S BELONGINGS HOLDER ON BASE

Usable area 295 in<sup>2</sup> (0.19 m<sup>2</sup>).



High-resistance materials and easy-to-disinfect finishes, designed to last and maintain high hygiene standards, compatible with products used in hospitals such as: hydrogen peroxide, quaternary ammonium, sodium hypochlorite and antiseptic alcohol.\*\*.

\*\*\* Recommended composition:

- Hydrogen peroxide 3% and 7.5% 200 ml of peroxide per liter of water.
- Quaternary ammonium 4 ml per liter of water.
- Sodium hypochlorite 4 spoonfuls for each liter of water approximately.
- Antiseptic alcohol 60-80% 7 parts of alcohol x 3 parts of water.





## MOBILITY AND EASE-OF-USE



### CENTRAL-LOCKING BRAKES

For secure and stable positioning, keeping the stretcher in place.



**Fixed adjustable IV pole** at head-end, with 2 service hooks.



### RETRACTABLE FIFTH WHEEL

Designed for enhanced maneuverability.






HORIZONTAL OXYGEN  
TANK HOLDER

**Minimum height 20.6 in**  
facilitates safe and easy patient entry and exit.

PATIENT SURFACE**	
Length	73 in (1854 mm)
Width	26 in (660 mm)

### MATTRESS (Optional with additional cost)



PAD COVER COLOR

TURQUOISE HOSPITAL  
EXTREME

Helps in the treatment and prevention of pressure injuries, high-density polyurethane (PU) foam core,  $26 \pm 2 \text{ Kg/m}^3$  (1,623 lb/ft<sup>3</sup>)\*.

### HOSPITAL EXTREME COVER

Made of nylon, permablok lacquer coating, light protection AATCC TM 16-3 / NTC 1479 / CFFA-2 a1 - Blue Wool, waterproof, anti-static ASTM D-257, bacteria resistance AATCC 147, JIS Z 2801, mildew resistant NO GROWTH (NG), invisible waterproofed seams, stitch resistant DIN 54301, flame retardant FMVSS 302 and latex free.

Ref. 740-570: 3 in MATTRESS FOR MAXCARE STRETCHER - HOSPITAL XTREME

Ref. 740-571: 4 in MATTRESS FOR MAXCARE STRETCHER - HOSPITAL XTREME

Ref. 740-572: 5 in MATTRESS FOR MAXCARE STRETCHER - HOSPITAL XTREME

STANDARD DESIGN MULTIPLE CONFIGURATIONS,  
WITH MOBILITY AND SAFETY IN MIND.

## MONITOR TRAY WITH CHART HOLDER\* Ref. 740-573

\*Optional Feature (Additional Charge).



## MOVEMENTS AND POSITIONS \*



Backrest



Horizontal




Trendelenburg



° Reverse Trend.

## GENERAL DIMENSIONS \*

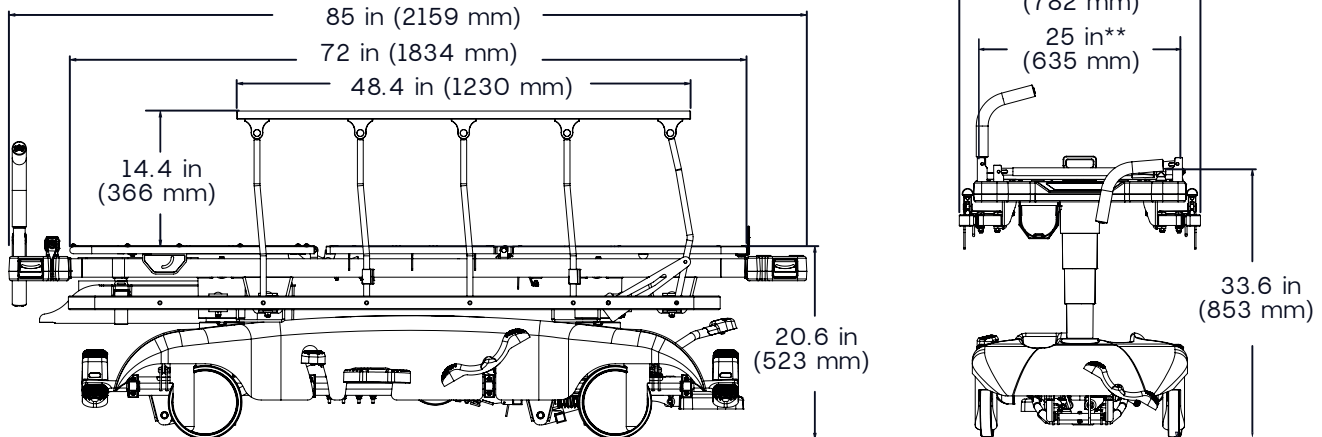

\* Movements performed by pneumatic and hydraulic mechanisms assisted by pedals on both sides of the stretcher.

\* Tolerances: degrees  $\pm 1^\circ$  dimensions  $\pm 0.5$  in (12.7 mm).

\*\*Platform width 25 in - Patient width 26 in.

## MANUFACTURING AND FINISHING

- Structure made of high-strength cold-rolled structural steel tubing.
- Finished in electrostatically applied epoxy-polyester powder paint resistant to oxidation with 3,000 hours of saline chamber.

## PACKAGING

- Container: 40HC - 39 units, 40STD - 26 units, 20STD - 12 units.



## STANDARD FEATURES

- Brake and steer pedals on 4 sides.
- Stationary knee section.
- Patient surface width 26 in.
- Retractable fifth wheel.
- Foldable push handles on head-end.
- 1 permanently mounted IV pole & 1 detachable telescopic IV pole, 2 stage with 2 hooks.
- Angle indicators, head of bed and trend/rev. trend.
- 8 in casters with covers.
- Pneumatic assisted backrest adjustment.
- Trendelenburg and reverse trendelenburg.
- 'Tuck away' collapsible safety side rails, with release lock.
- 4 integrated IV sockets.
- Drainage bag holders.
- 3 pedals to adjust height, trendelenburg and reverse.
- Central break system and fifth directional wheel.
- Total width of 30.7 in, facilitates entry and exit from hospital environments.
- Mattress retainer in foot section.
- Litter-mounted oxygen bottle holder.
- Supports for fastening straps in the back and feet sections.
- Integrated patient storage tray in base.

## SPECIFICATIONS

- 700 lbs. safe working load capacity.
- Patient surface width 26 in.
- Patient surface length 73 in.
- Backrest adjustment 0° – 90°.
- Height range 20.6 in – 33.6 in.
- Trendelenburg and reverse trendelenburg  $\pm 17^\circ$ .
- Total stretcher width 30.7 in.
- Total stretcher length 85 in.
- 3 hydraulic control pedals to adjust movements: height, trendelenburg and reverse trendelenburg on the sides of the stretcher and foot end.

## Optional Accesories (additional charge)

- Vertical oxygen tank holder. Ref. 730-675.
- Monitor tray with chart holder. Ref. 740-573.
- Detachable telescopic IV pole. Ref. 740-461.
- Fixed telescopic IV pole at the feet. Ref. 740-567-1.
- Foldable push handles at the feet. Ref. 740-566-1.
- 3 in. mattress with antibacterial and antistatic cover. Ref. 740-570.
- 4 in. mattress with antibacterial and antistatic cover. Ref. 740-571.
- 5 in. mattress with antibacterial and antistatic cover. Ref. 740-572.

## NOTES

- For accessories, consult your sales representative.
- Mubi Medical reserves the right to change specifications without prior notice.
- Warranty, consult your sales representative.



# TRANSFORMING SPACES, IMPROVING LIVES



Medellín – Colombia  
Calle 2 # 50-185  
(57) 604 285 2838

Cali – Colombia  
Calle 18 # 118-85  
(57) 602 489 6999



Date update 19/05/2025